## **Operation Instructions**

Light button Vibrating button 1) Power-on: Long press 0/0 for 2 seconds to power on or off. After boot-up, the device is default at: Temperature section 1 38°C.

o heat up at the targeted temperature depending on the environmental factors after power-on. It is recommended to use other functions (vibration and light) after power–on and heating. OF indicates the heating function is off. ③ Short press → : 4 kinds of light modes in a switch loop: near–infrared light (850nm) + red flashing light (630nm) →

near-infrared light (850nm) + yellow flashing light (580nm) →

near-infrared light (850nm). (Note: Near-infrared light is invisible to the naked eye);

④ Short press ★ : Vibration cycle adjustment: the mode is changed in sequence:

Section 1  $\$  Section 2  $\$  Section 3  $\$  Turn off the vibration,

### ⑤ Please adjust as appropriate functions to suit your personal needs. **Charging Instructions**

USB-C

charging port

display

USB-C

charging port

1. Before the first use, fully charge the device for about 2 hours; 2. When charging the product, a special USB cable is pr ovided to connect the TYPE-C socket of the device, the other part is connected with an arbitrary 5V adapter. 3. When the indicator light of the battery icon in the display

12. Do not use chemical solvents to clean the device, charging cable, etc.;

15. If the product has been used for a long time, the product should be adjusted to

The warranty provides free of charge repair within 1 year from the date of purchase. It also provides parts replacement including accessories within 6 months from the date of purchase. If there are any quality problems due to product defects,

our company shall be liable for the free repair, replacement of the product, and

Any damages due to accident, misuse, human error or transportation;

Any repairs/modifications made by any person other than authorized

> The product is misused which is not according the instructions manual;

> Any quality problems due to abnormal storage or protection conditions

> Any products purchased from unofficial or authorized consuming channels.

**Note:** The following disclaimer is not covered by this warranty.

Any losses is higher than the cost of the product;

**Environmental Protection and Recycling** 

➤ Unable to provide valid warranty card and registration;

The product has an internal installed rechargeable battery, which may

product to an officially designated battery recycling centre. You can also send

the discarded products back to our designated recycling department (please

contain sources of environmental pollution. Before disposing the product,

please be sure to turn off the power and take out the battery, or send the

contact the after-sales service centre for recycling information), and the

> Do not cause a short circuit of the battery in order to prevent a fire.

> Do not throw the battery into the fire in order to prevent explosion or

13. No disassembly. Otherwise, all consequences are at your own risk.

recommended to lower the temperature section or stop using;

the low (temperature) section.

replacement of the essential parts.

personnel of the manufacturer;

at the customer's workplace;

maintenance department will handle your case.

release of toxicants.

After-sales Service

One year warranty


14. If you feel that the temperature is too high / the skin turns red, it is

near-infrared light (850nm) + blue flashing light (460nm) →

the vibration icon will be off when the vibration is turned off;


panel keeps flashing, it indicates that it is charging; 4. When all indicator lights of the battery icon are on, it indicates that it is fully charged; 5. The product can be charged at any time and at any power level under particular environmental conditions;

6. When the battery icon shows that the battery is low, please charge it in time; 7. It is not recommended to fully use the battery until the battery level is low (The empty battery icon is on and flashing) to the product is shut down. 8. The product cannot be turned on when charging, please unplug the power supply before power-on.
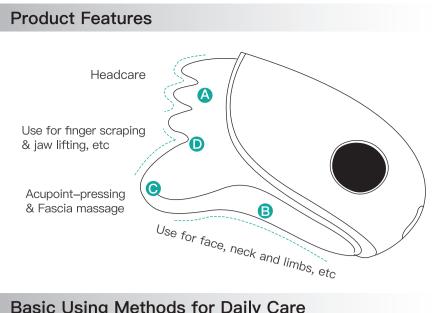
## **Basic Using Methods for Daily Care**
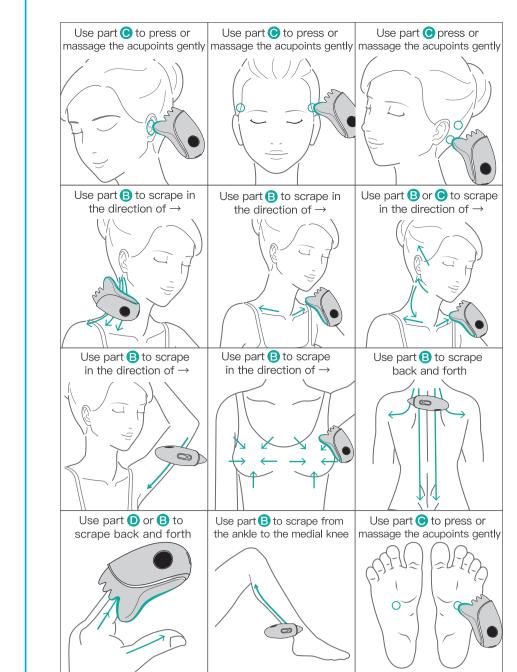
### Main Effects

√ To cooperate with acupoints for better effects

√ To enhance the absorption of skincare products into skin cells

- √ To relieve headache ✓ To make a wound heal faster √ To improve facial wrinkles and beautify skin √ To improve overall mental state √ To improve body microcirculation √ To be emotionally stable √ To improve metabolism and immunity ✓ To boost lymphatic drainage
- Applications of Energy Formula (Bianstone, Light Function)

Energy Type	Function Reference is made to "The Health Science of Far Infrared" and related optical literature.	Usage	Duration	Precautions
1. Bianstone releases far-infrared light (invisible light) after heating, negative ions, and ultrasonic waves are released when it vibrates.	Antioxidant, anti-inflammatory, improving microvascular circulation, pain relief, activating cells, and helping wound healing	1. Heating + Red light	20~40 mins	It is not recommended to use Heating + Blue Light at the same time
Near infrared     (invisible light)	Similar to far infrared	When the device is on, it continues working with every lights automatically.	20~40 mins	
3. Red light → ☆	Delay of aging, reducing pain and chronic inflammation (old wounds), improving skin wrinkles, removing freckles, enhancing hair growth, better sleep, improving mental state, and energizing.	Heating+Red light →	20~40 mins	1. It is not recommended for night time use. 2. Not suitable for people with epilepsy, purulent inflammation and tumor.
4. Yellow light - 次	Boosting lymphatic drainage, increase of adrenal secretion, lightening spots, improving self-confi- dence, concentration, memory, joy of emotion.	1. Heating+Yellow light	20~40 mins	Not suitable for people with hepatitis
5. Blue light 🔆	Relaxing nerves, reducing inflammation (new wounds), relief of palpitations, insomnia, bursting pain, oily skin control, reducing blackheads and acnes, and stabilizing emotions.	No heating, Blue light only	20~40 mins	Not suitable for people with hypothyroidism or depression.

### Maintenance, Cleaning & Storage

- 1. After use, please clean it with a damp towel, then dry it with a soft cloth; 2. Please store the product in moisture-proof, fire-proof, sun-proof, anti-pressure, anti-collision, anti-pollution, and dry and ventilated place;
- 3. Before and after use, please wipe the surface of the bianstones with rubbing alcohol; 4. After using the product, please turn off the power and place it in a safe place, and keep the product away to children;
- 5. When the product is unused for a long time, please place it in the storage bag. Do not expose the product under direct sunshine which is easy to make discolouration or deformation and damage. Meanwhile, please charge the battery once a month if the product has been unused for a long time, in order to avoid damage to the battery or dampness of the product; 6. Do not use this product for other purposes.t

### Specifications

Model: ER100

Name: Zikko Dr.Rock Energy Rejuvenator Main material: Natural red Bian Stone Auxiliary heating: Graphene film Temperature adjustment: Section 1 38°C, Section 2 45°C, Section 3 50°C (Note: A deviation of ±5°C will be occurred depending on environmental and usage conditions)

Vibration: 3-speed adjustment **Light wavelength**: Far infrared 6 ~ 14 μm, near infrared 850nm, red light 630nm, yellow light 580nm, blue light 460nm Visible light flicker frequency: 7.83 Hz Ultrasonic frequency: About 650 KHz

Power input: USB-C 5V/1A Lithium battery: 3.7V/1000mAh, 3.6Wh Battery life: 1 hour (it may vary according to usage) Charging atime: About 2 hours Size: L 11\* W 7.6\* H 3.7 CM

## Package & Accessories

 Energy rejuvenator x1 Product storage bag x1

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Weight: 165g

Color: White

 User manual x1 USB charging cable x1

### Abnormal Fault & Trouble-clearing

LED polarized light

Fault Description	Reasons for Abnormal Fault	Trouble-clearing
The product is not turned on.	Please check the battery level.	Please charge the battery. (Please charge at least 2 hours if it has been unused for a long period of time.)
Th product shutdown when using	The battery is insuffient.	Please charge the battery.
The product is failed to be charged or has no charging display.	Please check the connection with the power supply and plug it fully.     Please check the battery level of the power supply.	Please re–plug the power or change a new adapter.     Please check the connection with the power supply.
No vibration response	The vibration part is not connected with the power supply or is damaged.	Please contact after sales service centre for specific reasons.
No heat	Please wait for 5 ~ 15 mins. / Check the connection of heating	Please contact after sales service

### centre for specific reasons. films with power supply. Note: For other abnormal issues, please contact customer service centre immediately.

### Important Safety Instructions

- 1. Please read the instruction manual carefully before using the device; 2. This product is not a medical device, it is not recommended for medical use; 3. The product has a heating surface and is not suitable for people who are not sensitive to heat, infants or people with reduced physical, sensory or mental capabilities; 4. It
- is not suitable for children to use, and it requires parents or guardians to set it up in advance, or children have been given supervision or instructions; 5. When using, please pay attention to the temperature of the heated bianstones, choose a suitable temperature section, and choose the appropriate temperature
- which the skin can be generally accepted. Do not overheat in order to prevent burns; 6. When the product has odor residue after using with essential oil, please wipe with a clean damp cloth in order to keep the product clean. It is recommended to use 75% alcohol for disinfection before use in order to avoid infection; 7. The bianstone part of this product is made of natural xuan-huang bianstone, the
- appearance of fine texture is a normal phenomenon, which has no effect on the use of bianstone; 8. Please do not use the product during charging;

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9. Please shut down before unuse; 10. This product cannot be dropped or hit from high altitude, it has to be handled with care, and it is fragile;

### 11. The device is not waterproof, do not use it with wet hands and feet or in a Toxicant name and content in the product humid environment;

	Toxicants or elements					
Component Name	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent chromium (Cr(VI))	PBBs (PBB)	PBDEs (PBDE)
Plastic Case	0	0	0	0	0	0
Bianstones	0	0	0	0	0	0
Rubber parts	0	0	0	0	0	0
Circuit board components	0	0	0	0	0	0
Cables	0	0	0	0	0	0
Motor	0	0	0	0	0	0

This form is compiled in accordance with SJ/T 11364. O indicates that the content of toxicants in all materials of the part is

below the limited requirement specified in GB/T 26572. X indicates that the content of toxicants in at least one of a certain materials of the part exceeds the limited requirement specified in the GB/T26572.

1. During the environmental protection period of usage, please use the product in accordance with the instruction manual. 2. The environmental protection period of usage does not cover consumable items, batteries, motor.

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